

Product Specifications

Device	Weight	Height	Width	Distance from Pylon	Depth	Pump Space Volume
iQ Vacuum Device	271 g	4.06"/ 103mm	2.34"/ 59mm	2.83"/ 71.8mm	3.89"/98.9mm	19.2 in ³

Part Description	Part Number
iQ Vacuum KIT	VAC-00-A1001-KT

Part Description	Part Number	Qty	Included in Kit	Available Separately?	Warranty
iQ Vacuum System	VAC-00-A1001-00	1	Yes	No	24 months
Pylon Bumper, Foam Tape, 4"	OTS-32-00008-00	1	Yes	Yes	N/A
Tube, 1/16" ID X 12"	OTS-39-00030-00	12 in	Yes	Yes	N/A
AC/DC Wall Mount Adapter	OTS-62-00005-00	1	Yes	Yes	6 months
Fitting, 1/8" 90 Degree Elbow Barb	OTS-69-00010-00	1	Yes	Yes	N/A
Filter Assembly X 4	VAC-00-A1011-00	4	Yes	Yes	12 months
1.142" Clamp Bushing	VAC-00-P0101-01	1	Yes	Yes	N/A
1.282" Clamp Bushing	VAC-00-P0101-02	1	Yes	Yes	N/A
1.370" Clamp Bushing	VAC-00-P0101-03	1	Yes	Yes	N/A
iQ Vacuum User Guide, Clinician	CG-03-007	1	Yes	Yes	N/A
iQ Vacuum User Guide, Patient	CG-03-008	1	Yes	Yes	N/A
iQ Vacuum IFU	IFU-03-004	1	Yes	Yes	N/A

Focused Solutions

Freedom Innovations is solely focused on developing world-class lower-limb solutions in close collaboration with prosthetists and people with amputations. These solutions encompass technological innovations, service, training, educational resources and consultative support to help ensure user satisfaction.

www.freedom-innovations.com

Global Headquarters

3 Morgan
Irvine, CA 92618
toll free: 888-818-6777
phone: 949-672-0032
fax: 949-672-0084
www.freedom-innovations.com

European Head Office

Jaargetijdenweg 4
7532 SX Enschede
The Netherlands
toll free: 00800 - 2806 2806
phone: 31 (0) 53 20 30 300
fax: 31 (0) 53 20 30 305
www.freedom-innovations.eu

USA Manufacturing & Returns

425 East 400 North
Gunnison, UT 84634




#freedominnovations



Intelligent Design
for Precision Performance

Candace,
Southern
California



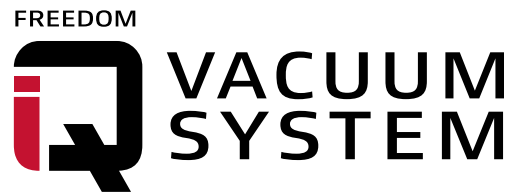

Intelligent Design for Precision Performance


For prosthetists who seek to deliver new gains in clinical outcomes, efficiency, control and predictable performance, the next generation of electric vacuum has arrived.



The No Compromise Vacuum Solution

iQ is the first electric vacuum that overcomes the inherent limitations of traditional vacuum while providing the combined benefits of electric vacuum in an all-in-one device, offering:

- Automatic adjustment
- Unparalleled quietness
- User and prosthetist apps for streamlined setup, usage and cadence reporting
- External Set-up design for trial and retrofitting option
- Customizable settings
- IP67 rating

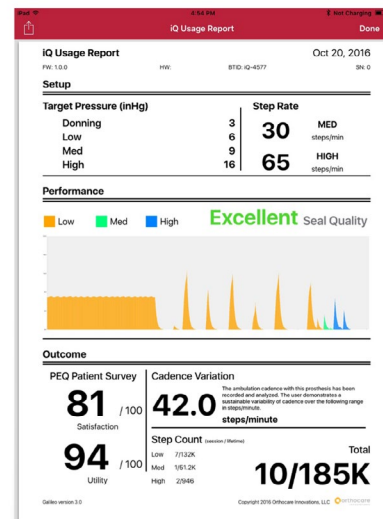

Microprocessor-Controlled system

Microprocessor-Controlled System

Real-time responding sensors detect user's changes in activity (walking vs. standing vs. sitting) to automatically adjust the vacuum pressure. Built-in leak detection capability provides reliable consistent socket comfort.

IP67 Rating

iQ Vacuum is protected against ingress of dust and submersion in water at a depth of 1 meter for up to 30 mins, allowing user to withstand occasional wet environments.



Unparalleled Quietness

iQ operates at the most quiet noise level available in a normal environment, eliminating user embarrassment and annoyance.

External Setup Design and Build Height Flexibility

iQ requires minimal space on standard Pylon to allow for flexibility with build height. The outside socket design allows for trailing the product as well as retrofitting options.

User and Prosthetist App

Streamlined set up, monitoring, customization and cadence report out are a touch away with the Clinician and Patient App.